

# ULTRASOUND MAGIC 6000

*High Resolution Images*

*Large  
Screen*

*Sharp  
Images*




## Ultrasound

# HIGH RESOLUTION IMAGING MAGIC 6000

---

*Partnering with veterinarians in clinical practice has helped develop the products that clinicians need. The MAGIC 6000 provides crisp and clear diagnostic images for veterinarians.*

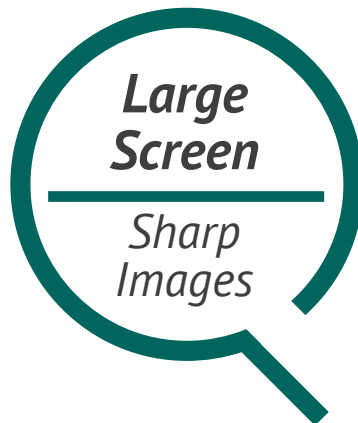

### Clear Image

Pleasant image presentation even at a longer viewing distance.

### Simple Archive System

iStation™ helps you find the correct patient data quickly and efficiently.

### Easy to Use

The combined keyboard and menu navigation controls make it uncomplicated for beginners to find their way.

### Ideal for Small Animal Practice

The second transducer connector allows the use of multiple probes during the scan.

---

*The MAGIC 6000 is a reliable companion in my mixed practice. With it, I can examine small and large patients easily.*

---

*compact*



## ULTRASOUND MAGIC 6000

- ▶ Anatomic specific measurements: Gastrointestinal and reproductive organs
- ▶ Standard measurements: distance, circumference, diameter, volume
- ▶ Imaging modes: B, B+B, M, B+M
- ▶ Transducer ports: 2
- ▶ Large internal memory (320 GB)
- ▶ iStation™ patient management, image and data storage
- ▶ Cine Mode with storage function (5 seconds)



- ▶ Screen size: 15" LCD
- ▶ iBeam™: Reduces random noise and enhances border differentiation
- ▶ iClear™: Speckle reduction imaging
- ▶ iZoom™: Automatically maximises picture to full screen
- ▶ iTouch™: Automatic gain compensation
- ▶ 8-layer TGC: Balances the signal strength of different layers
- ▶ Image export in Window format
- ▶ Raw image data storage
- ▶ 4 USB-ports
- ▶ Dimensions (in cm): W 41.5 x H 39.6 x D 47.6
- ▶ Weight: 8 kg (without battery)

**303720**

## Ultrasound

# ACCESSORIES FOR MAGIC 6000

*The technology behind these transducers technology provide superior image quality. Transducers are available in several versions, so we can offer a solution that meets the needs of your clinical interest.*



### TRANSDUCER 6.5 MHZ MICRO-CONVEX

- Frequency: 4.0/6.0/7.5/8.5/  
H8.0/H9.0 MHz

**303710**



### TRANSDUCER 7.5 MHZ LINEAR T-SHAPED

- Frequency: 5.0/7.5/8.5/10.0/  
H8.0/H10.0 MHz

**303711**



### TRANSDUCER 5.0 MHZ LINEAR RECTAL

- Frequency: 4.0/5.0/6.0/7.0/  
H6.0/H7.0 MHz

**304021**



### TRANSDUCER 7.5 MHZ LINEAR RECTAL

- Frequency: 5.0/7.5/8.5/10.0/  
H8.0/H10.0 MHz

**304022**



### CARRY CASE FOR MAGIC 6000

- To help protect during transport


**303725**



### TROLLEY FOR MAGIC 6000

- Helps transport the machine  
between rooms

**303723**

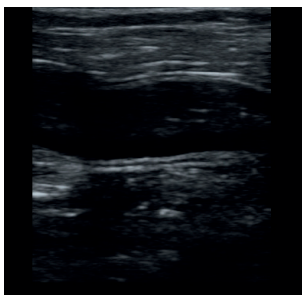

*Duodenum of a small dog with  
a linear transducer*



*Adrenal gland of a small dog with  
a convex transducer*

### BATTERY FOR MAGIC 6000

- 2 hour battery life

**303721**