

Life's better with **HYDRO PHYSIO**[®]
UNDERWATER TREADMILLS




CANINE TREADMILLS AND POOLS



Delivering Canine Care
with Confidence




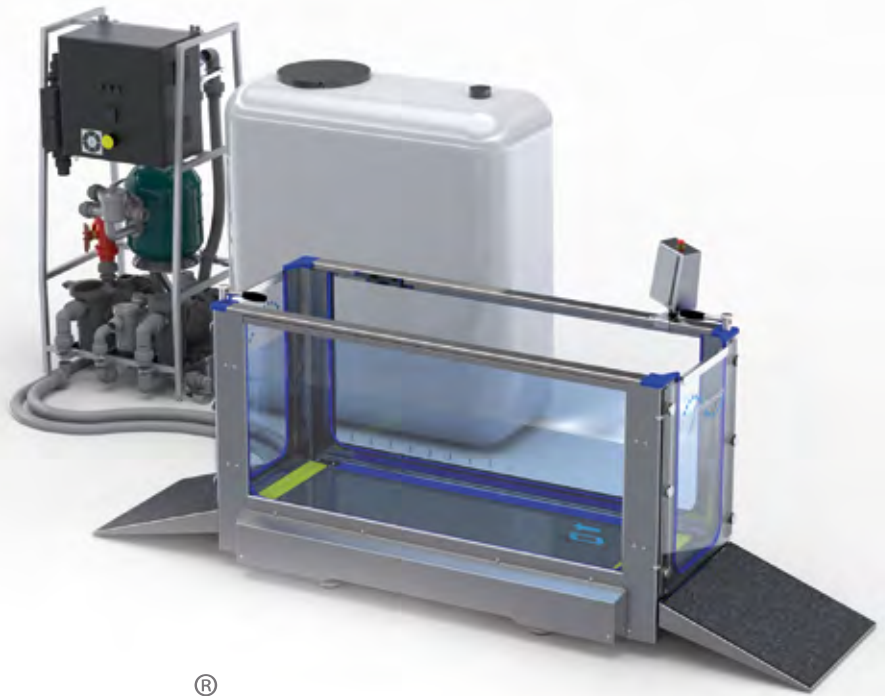

HYDRO PHYSIO[®]

Hydro Physio Underwater Treadmills and Pools - Delivering Canine Care with Confidence

With more than 20 years of designing and manufacturing canine underwater treadmills and pools, we support Hydrotherapists in delivering enhanced health, recovery and fitness outcomes for their patients.

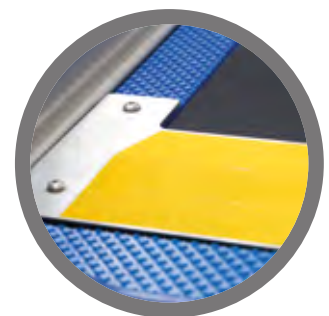



HYDRO PHYSIO®

Hydro Physio leading the way

Based in Shropshire, England, Hydro Physio have been supplying quality underwater treadmills and pools for over 20 years. In our recently expanded 200,000 Square Foot facility, the design, manufacturing, sales and service teams are all under one roof, facilitating unrivalled customer support.

Innovative design ensures Hydro Physio treadmills have cutting edge features, whilst ISO compliant manufacturing assures a high quality product. The experienced service team offer expert guidance via telephone support, and fully trained service engineers will visit your facility when required.

Over 1,000 Hydro Physio systems have been installed in over 50 countries worldwide.



Thank you to the team at Hydro Physio for their excellent service from ordering to installation. Everyone was friendly and helpful in explaining the process and setting up the treadmill. My spaniel, Looney, very much enjoyed having the guys around to install the treadmill. Thanks again, I am very pleased with the final product!



Ellen Winter
Veterinary Physiotherapist

Discover Hydro Physio Canine Treadmill

Designed with input from expert hydrotherapists and veterinary professionals, our treadmills are manufactured to provide patients with water-based rehabilitation.

Our treadmills come with a **two-year Parts and Labour warranty***, giving you peace of mind with your investment.

●●●● Key benefits of our treadmill include

● Exceptional build quality

- Manufactured from medical-grade stainless steel
- 360-degree glass for all-round patient visibility

● Practitioner-friendly features

- Side platforms for easy access
- Remote control from portable devices
- Remote service diagnostics
- Adjustable control panel on chamber
- Powered incline for changing workout intensity

● Excellent safety and accessibility

- Walk-through entry/exit door with a patented safety latch
- Anti-slip surfaces in the chamber and on the ramps
- Reverse belt direction for flexibility (max speed applied)

●●●● Additional options for practitioner convenience

- Side deck
- Therapy stool
- Isolation plate

●●●● Multiple installation options

- Stand-alone: Free-standing
- Stand-alone: In-ground
- Pool-fed

Arrange a site survey

Our expert team can visit your facility to identify your specific installation requirements. Following your free site survey, we will provide you with detailed layout drawings.



Discover Hydro Physio Pools

If you are looking for a cost-effective swimming pool that facilitates rejuvenating underwater workouts, our free-standing, above-ground pools offer the perfect solution.

Designed for easy installation and relocation, **Hydro Physio pools** are ideal for therapists looking to provide patients with a convenient, engaging exercise option. Available in multiple sizes from 18 m³ - 56 m³, they are suitable for a wide selection of facilities.

●●●● Key benefits of our free-standing pool include:

● **Made from high-tenacity polyester fabric**

● **PVC coating ensures safety and cleanliness**

● **Easy to assemble and dismantle**

● **Pool frames are constructed from robust and rigid anodised aluminium**

* UK based customers only (first years service needs to be completed to maintain the additional 2nd year warranty).

The Sectors We Serve

Veterinary Community

Provide patients with joint-friendly exercise that improves recovery, pain management and overall quality of life.

Hydrotherapy

Enhance pain management, mobility and rehabilitation with our underwater treadmills and pools.

Animal Welfare

Help vulnerable animals regain strength and mobility in a safe, controlled environment.

MOD and Police

Support working dogs in maintaining exceptional levels of fitness and agility.



Treadmill Specification

Maximum load	90kg (198lbs) without buoyancy
Maximum water depth	0 - 600mm (24")
Speed range	0.2 - 12kph (0.12 - 7.5mph)
Fill / empty times	Approx 2 minutes to max depth
Temperature range	From ambient - 40°C (104°F)
Treadmill dimensions (Excluding display) E	
Excluding ramps	L 2150mm (85") x W 860mm (34") x H 1265mm (50")
Including ramps	L 3435mm (135") x W 860mm (34") x H 1265mm (50")
Complete system weight	
Dry	650kg (1433 lbs)
Full	2150kg (4740lbs)
Water tank capacity	1500 litres (330 gals, 396 US gals)
Max distance between treadmill and plant equipment	8m (314")
Power	240V single phase 50hz 16A (60hz option available)


Pools Specification

MODEL	EXTERNAL MEASURE MM	INTERNAL MEASURE MM	APPROX VOLUME
MAX-1	5400 x 3650 x 1200	4750 x 3000 x 1200	18m ³
MAX-2	6650 x 3650 x 1200	6000 x 3000 x 1200	23m ³
MAX-3	9150 x 3650 x 1200	8500 x 3000 x 1200	32m ³
MAX-4	10400 x 4900 x 1200	9750 x 4250 x 1200	50m ³
MAX-5	11650 x 4900 x 1200	11100 x 4250 x 1200	56m ³





Life's better with **HYDRO PHYSIO**[®]
 UNDERWATER TREADMILLS

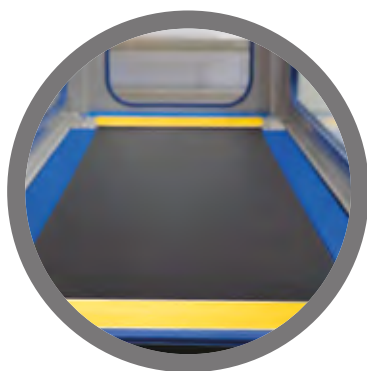

●●●● **Get in Touch**

+44 (0) 1952 885112

canine@hydrophysio.com

hydrophysio.com



OHS 672496



FM 25162

Delivering Canine Care
 with Confidence

